


GREEN LASER PHOTOCOAGULATOR



MERILLAS 5320α

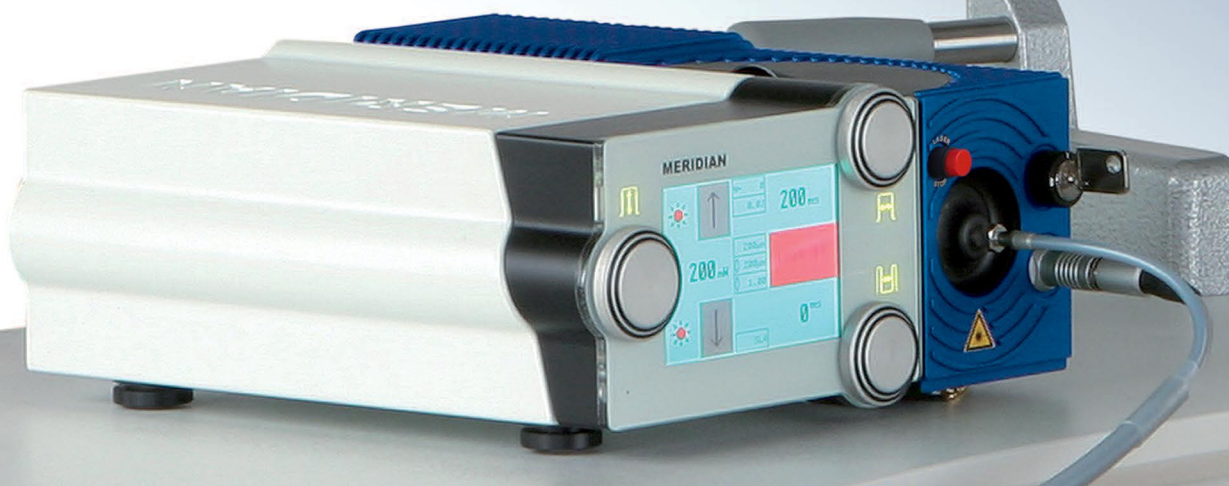

MERIDIAN
1937



MERILAS 532 α

GREEN LASER PHOTOCOAGULATOR

THE PROFESSIONAL FREQUENCY DOUBLED SOLID STATE GREEN LASER PHOTOCOAGULATOR FOR PHOTOCOAGULATION WITH SLIT LAMP ADAPTER, OPHTHALMOSCOPE, AND ENDOPROBE.

The compact MERILAS 532 α for reliable and precise photocoagulation procedures. The unique thermo electrical cooling system «Smart Cool» ensures absolutely quiet operation, without the irritating sound of a cooling fan.

The detachable user interface with the touch screen glass technology, ensures the highest convenience in operating and cleaning. The «Fast Start» software allows a quick startup time.

The α line of ophthalmic lasers were designed for optimised ease of use, treatment efficiency, and safety.

With its compact size and weight it is the ideal unit for the environment where high performance combined with low space requirements is mandatory.

Even with its compact size and footprint the MERILAS 532 α is equipped with a detachable user interface with colour glass touch screen and graphic user guidance.

Because of its ergonomic design the MERILAS 532 α has been nominated for the Swiss Design Prize 2007.



The easy to operate detachable user interface can be placed up to 1 meter away from the laser

MERILAS 532 α – the green 532 nm laser for

- All retinal photocoagulation procedures
- Iridotomy
- Laser Trabeculoplasty

MERILAS 532 α – the smart choice

- Small, portable, and light weight unit
- Ease of use due to graphical user interaction
- «Fast Start» software for short time
- Quiet operation without the noise of a cooling fan
- User interface detachable for user convenience with glass touch screen, easy to clean and highly resistant
- Low maintenance requirements and long service life due to diode pump technology



Ergonomic, plenty of space available and no interfering objects between the slit lamp front lens and the patient's eye.

Slit lamp adapter

- Best suited for performing retinal photocoagulation procedures, iridotomy and Laser Trabeculoplasty
- Directly linked to the illumination path of the HAAG-STREIT BM 900 and BQ 900 slit lamp – no additional obstacle between slit lamp front lens and patient's eye
- Safe treatment due to internal beam path – no external parts which can be misaligned

Laser Indirect Ophthalmoscope LIO 500

- Ideal for treating patients in a supine position
- Fully enclosed optic for optimal protection



The slit lamp adapter mounted directly in the illumination path of the slit lamp...

Endoprobes

- For endophotocoagulation procedures
- Straight and angled probes available
- Safety filters for full protection available

User protection filters

- Unique «Clear View» user protection filter for highest visibility and low colour distortion available for HAAG-STREIT slit lamps and a large range of surgical microscopes
- Active (motor driven) user protection filter available for HAAG-STREIT slit lamps and a range of surgical microscopes.

Optional Accessories

- Mounting set for MICRORUPTOR 6
- Safety goggles
- HAAG-STREIT laser contact lenses



...and mounted in the slit lamp adapter holder.

MERIDIAN has been designing, manufacturing and maintaining ophthalmic lasers for more than twenty-five years. Since the introduction of the first Nd:YAG laser by Prof. Fankhauser, University of Berne Eye Clinic in 1982, Thun has become the place where superior ophthalmic lasers are designed and built. A worldwide distributor network ensures comprehensive service support throughout the world.

See www.meridian.ch for details.

MERIDIAN IS PARTNER OF



MERIDIAN

1937

SPECIFICATIONS

MERILAS 532 α

Part No. 100491

Treatment Laser Wavelength: 532nm

Treatment Laser Class: IV

Supply Voltage: 90 – 260 VAC 50/60 Hz

Protection Class: I

Power Consumption: 250W

Cooling: Thermo electrical

Maximum Output Laser Power: 2000mW

Pulse Duration: 10ms – 5s

Delivery Systems: Slit lamp adapter,

Ophthalmoscope LIO, Endo fibre (SMA)

Aiming Beam Laser Wavelength: 635 nm

Aiming Beam Laser Power: 0.01 – 1mW adjust.

Aiming Beam Laser Class: II

Dimensions: 105 (H) x 220 (W) x 250 (L)mm

Weight: 7 kg

Packaging: Protective carrying/storage case (incl.)

OPTIONAL ACCESSORIES

SLIT LAMP ADAPTER for HAAG-STREIT

Part No. 340835

Suitable for HAAG-STREIT BQ and BM 900 slit lamps, directly linked into the illumination path of the slit lamp, spot size: 50 μ m to 5000 μ m

USER PROTECTION FILTERS Part No.

Passive for BQ 900 340843

Active for BQ 900 340850

Passive for BM 900 340845

Passive for Zeiss/Moeller Wedel 340833

Active for Zeiss/Moeller Wedel 340849

Passive for Leica 340846

Active for Leica 340848

SLIT LAMP ADAPTER HOLDER

Part No. 311070

SLIT LAMP ADAPTER UNIVERSAL

Part No. 340847

Suitable for mounting on the tonometer pin of slit lamps, spot size: 50 μ m to 1000 μ m

LASER INDIRECT OPHTHALMOSCOPE

LIO 500

Part No. 340980

Including illumination-kit, spot size: 1 mm (without laser lens)

LASER LENSES HAAG-STREIT Part No.

CGIL, Iridectomy 7200194

CGAL, Gonioscopy 7200195

CGPL, Capsulotomy 7200196

CGRL, Retina 7200197

CGVL, Vitrectomy 7200198

OTHER HAAG-STREIT ACCESSORIES

Available upon enquiry, MERIDIAN is partner of HAAG-STREIT INTERNATIONAL



The Merilas 532 α in combination with the MICRORUPTOR 6 makes a combination unit with the ergonomical Nd:YAG laser available.



Caution: Federal law restricts this device to sale by or on the order of a physician.

Your distributor:



MERIDIAN

Biergutstrasse 7

CH-3608 Thun, Switzerland

Phone +41 (0) 3333 4 1111

Fax +41 (0) 3333 4 1119

info@meridian.ch, www.meridian.ch

MERIDIAN

1937