

Totaalleverancier van lasers en medische technologieën



méér dan lasers ...



: HEINE BETA® LED F.O. Otoscopes

The best and brightest HEINE Otoscopes based on 70 years of experience and exclusive LED^{HQ} Technology



[HEINE BETA® 400 LED F.O. OTOSCOPE]



HEINE BETA® LED F.O. Otoscopes: Sharpest images, perfect colour rendering and brightness

Exclusive 4.2x or 3x magnification. More detail of anatomical structures and easy recognition of foreign bodies.

Fibre Optic Illumination. Ensures homogeneous, very bright illumination and an unobstructed view of the ear canal and tympanum.

Exclusive continuous brightness control between 100% and 3% (patent pending). **Optimal light intensity can be individually adjusted for every examination situation.**

The permanently attached optical system which can be flipped open/swivelled out of the way enables instrumentation and easy cleaning of the instrument head.

With an all metal housing and a scratch resistant glass window, BETA LED F.O. Otoscopes **guarantee a lifetime of maintenance-free use** and high-quality patient assessments.

5 Year Guarantee. HEINE makes no compromises in manufacturing high quality medical instruments. Our commitment to vertical integration in manufacturing means we control all aspects of our instrument quality.





HEINE LEDHQ: Light that shows colours as they are



LED light sources impress with their long lifetime and energy efficiency. There is no need to ever exchange the LED. The quality of LED illumination can however vary greatly compared to XHL illumination.

With the exclusive HEINE LEDHQ technology you can enjoy all benefits of LED technology, with no compromise in illumination quality. This is achieved by a high Colour Rendering Index (CRI) of 97 (red is red, blue is blue).

BRIGHT AND HOMOGENEOUS ILLUMINATION

with typ. 77,000 lux (with 5 mm tip).

HIGH COLOUR RENDERING INDEX

Colour temperature: typ. 3,500 K, Colour Rendering Index (CRI) > typically 97, special index for red colours (R9) > 93.

CONSISTENT LIGHT OUTPUT AND PERFORMANCE

throughout the whole working life of typ. 100,000 hours.

EFFICIENCY

with the BETA4 rechargeable handle a longer examination time of approx. 3x.

[HEINE BETA® 200 LED F.O. OTOSCOPE]



[ORDERING INFORMATION]

BETA400 LED F.O. Otoscope Set with BETA4 USB rechargeable handle

The set includes the BETA400 LED F.O. Otoscope; 1 set (4 pcs.) reusable tips and 5 of both 2.5 and 4mm Ø AllSpec disposable tips, **BETA4 USB rechargeable handle** with **USB cord and plug-in power supply,** hard case, handle rest


B-143.28.388



BETA 400 LED F.O. Otoscope Set with BETA 4 NT rechargeable handle

The set includes the BETA400 LED F.O. Otoscope; 1 set (4 pcs.) reusable tips and 5 of both 2.5 and 4 mm Ø AllSpec disposable tips, **BETA4 NT rechargeable handle with NT4 table charger**, hard case

B-143.24.420

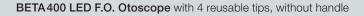

B-008.11.401



BETA 200 LED F.O. Otoscope Set with BETA 4 USB rechargeable handle

The set includes the BETA200 LED F.O. Otoscope; 1 set (4 pcs.) reusable tips and 5 of both 2.5 and 4 mm Ø AllSpec disposable tips, **BETA4 USB rechargeable handle** with **USB cord and plug-in power supply,** hard case, handle rest

B-141.28.388



BETA 200 LED F.O. Otoscope Set with BETA 4 NT rechargeable handle

The set includes the BETA200 LED F.O. Otoscope; 1 set (4 pcs.) reusable tips and 5 of both 2.5 and 4 mm \emptyset AllSpec disposable tips, **BETA4 NT rechargeable handle with NT4 table charger,** hard case

B-141.24.420



B-008.11.501





For additional product configurations and further HEINE Quality Products please visit **www.heine.com** or give us a call!

We reserve the right to change specification without notice.

© 2017 HEINE Optotechnik GmbH & Co. KG. All rights reserved

Your distributor

07/17. A-000.00.266 e













De Grift 20 7711 EJ, Nieuwleusen