

Totaalleverancier van lasers en medische technologieën



méér dan lasers ...



This is how confidence feels



LightSheer® QUATTRO™

Feel the confidence that comes with Lumenis Laser Hair Removal




We don't just build technology, we build confidence.

What we do changes the way people feel about themselves.

We don't just remove unwanted hair, we remove doubts and worries and the barrier between people and what they can be.

We create, design and manufacture confidence.

This is how **confidence** feels

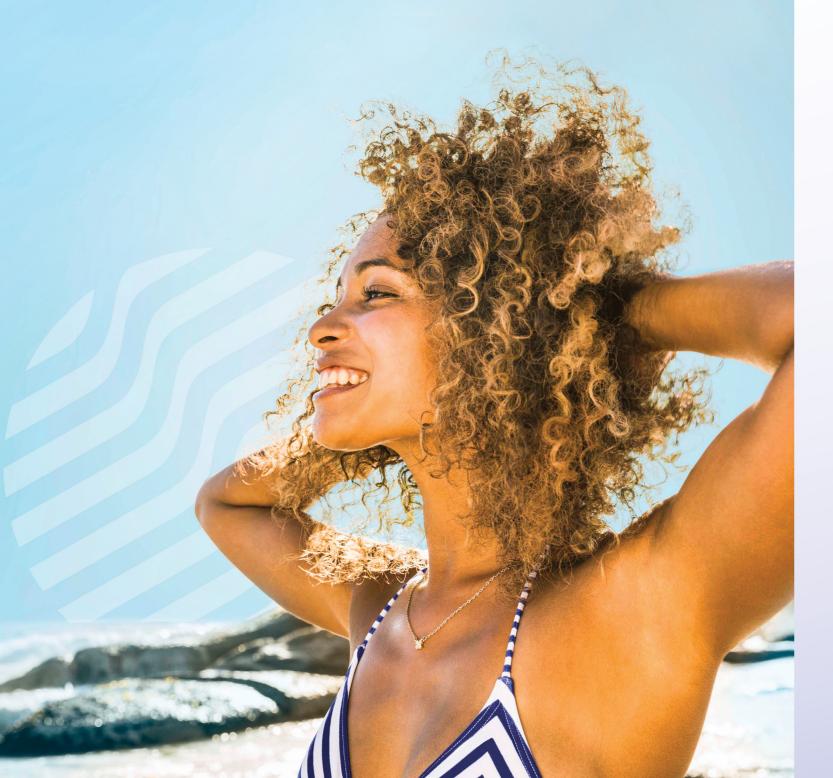

Satisfying patients, growing business.

Laser hair removal is one of the most sought after non-invasive aesthetic treatments for women and men alike. More and more people from different ethnic backgrounds seek permanent hair removal, they are looking for a treatment that is effective, fast and comfortable.

Laser hair removal draws patients to your clinic and serves as an excellent stepping-stone for them into the world of Aesthetic Dermatology.

With LightSheer QUATTRO you can customize patient treatments according to each individual's hair type, skin type (including tanned and dark skin) and lifestyle - and therefore access a larger demographic.



Any skin type with 805nm & 1060nm wavelengths



Upgradeable and scalable



Portable with a compact footprint



Treatment versatility pigmented & vascular lesions



You can be **confident** in our **technology**.

Great clinical results begin with the LightSheer QUATTRO diode laser empowered by two advanced technologies: High Speed vacuum assisted technology and ChillTip™ integrated contact cooling technology

Treat ALL SKIN TYPE with either 805nm or 1060nm wavelengths

High Speed Vacuum Assisted Technology



Efficacy²

Laser energy delivered more efficiently.



Safety²

Lower fluence is required per pulse.



Decreased pain¹

Unique mechanism reduces pain during treatment.



Treatment time^{1,2}

Rapid coverage rate without gel or anesthesia.

How it works







- 1 Ibrahimi OA, Kilmer SL. Long-term clinical evaluation of an 800-nm long-pulsed diode laser with a large spot size and vacuum-assisted suction for hair removal. Dermatol Surg. 2012
- 2 Xia Y, Moore R, Cho S, Ross EV. Evaluation of the vacuum-assisted handpiece compared with the sapphire cooled handpiece of the 800-nm diode laser system for the use of hair removal and reduction. J Cosmet Laser Ther. 2010; 12: 264-268.
- 3 Campos, V. B., Dierickx, C. C., Farinelli, W. A., Lin, T. Y. D., Manuskiatti, W., & Anderson, R. R. (2000). Hair removal with an 800-nm pulsed diode laser. Journal of the American Academy of Dermatology, 43(3), 442–447.

ChillTip Integrated Contact cooling Technology



High Speed handpiece

Offers a unique pain reduction

quickly and comfortably.

mechanism and works without gel.

Ideal for treating larger body areas,

22x35mm

Efficacy^{2,3}

Compression of the skin reduces competing chromophores.



Safety³

Epidermal protection through continuous sapphire contact cooling.



Patient comfort³

Continuously cooling significantly increases comfort.



XC

12x12mm

Providing a solution for large and small areas requiring higher fluence levels.









E

9x9mm

Designed to treat small areas that require maximum precision.

4 LightSheer QUATTRO LightSheer QUATTRO

We're confident in our results, so you can be too.

The LightSheer technology by Lumenis has been validated in numerous clinical studies and peer-reviewed articles, and our customers the world over have continuously reported high levels of patient satisfaction and repeatable results.

What sets LightSheer®

apart is the confidence it

gives providers in treating

patients comfortably and

with a high degree of

efficacy. It's an overall

provider and the patient.

Daniel P Friedmann, M.D.

win-win for both the

Laser hair removal, 805nm wavelength

Before



Pigmented lesions, 805nm wavelength





Wrinkle reduction, 1060nm wavelength





Laser hair removal, 1060nm wavelength





Photo courtesy of: 1. Mitchel Goldman, 2, 3. Lumenis Ltd. M.D. 4. Sushil Tahiliani, M.D.

Warnings and Risks:

Treatment should not be attempted on patients with the following conditions in the treatment area: any active infections, inflammatory skin conditions, Dysplastic nevi, tattoos, active cold sores, open lacerations or abrasions, viral, fungal, or bacterial diseases or a history of post inflammatory hyperpigmentation. The most common immediate responses of laser treatment for hair removal, ightening/darkening of pigment, vascular lesions, and wrinkles pigmented lesions are erythema, edema, perifollicular edema and perivascular edema, vascular bleaching, hyperpigmentation and hypopigmentation or textural changes of the treated areas.

The complete list of contraindications and risks can be found in the LightSheer QUATTRO user manual.

Your journey with Lumenis has just begun.

We want your business to grow with LightSheer QUATTRO. We provide a support package which ensures that you have all the training you need and recieve the most up to date clinical materials. Our direct to consumer marketing activities support the brand, keeping LightSheer the treatment of choice when it comes to laser hair removal.





Thousands of patients search for treatment providers every month. Help them find your clinic by registering your practice on the Lumenis Practice Locator. We run campaigns through influencers, social media and viral videos to promote this directory.



Lumenis online

Consumers have many questions about aesthetic. Our patient-facing websites will give them the re-assurance and confidence they need to move forward with the treatments.



Team Training

We install the system and train your staff, so that you are able to assess skin types and conditions, creating the perfect treatment plan for each individual patient.

Our expert clinical training means they will handle consultations and treatments with confidence.



Your Most Valuable Resource

Your dedicated representative is an email, text, or phone call away. We want you to feel part of the Lumenis family, so please feel free to contact us any time. We'll always do our best to get back to you within 24 hours.



Brand Your Clinic

As a Lumenis client you will have access to the Lumenis Physician Portal. Here, you can find materials such as brochures and posters. You will also have access to a wealth of clinical and product information.



Extensive Service Network

We are committed to providing high quality, comprehensive service for medical equipment repair and maintenance. Lumenis service provides you with peace of mind, so you can focus on patient care, rather than on the equipment.

6 LightSheer QUATTRO LightSheer QUATTRO 7

Specifications

LightSheer QUATTRO	HS 22 X 35 mm		XC 12 X 12 mm		ET 9 X 9 mm	
Laser Type	AlGaAs Diode Arrays		AlGaAs Diode Arrays		AlGaAs Diode Arrays	
Wavelength	805 nm	1060 nm	805 nm	1060 nm	805 nm	1060 nm
Repetition Rate	Up to 3 Hz		Up to 3 Hz		Up to 3 Hz	
Pulse repetition	1-3 pulses		1 Pulse		1 Pulse	
Fluence	805 nm: 4.5-12 J/cm ²	1060 nm: 4.5-14 J/cm ²	10-40 J/cm ²		10-100 J/cm ²	
Pulse Width	30-400 ms		5-400 ms		5-400 ms	
Peak Power	2310 W		2200 W		1620 W	
Preset Parameters	Yes		Yes		Yes	
Comfort Management	Vacuum assisted (up to 18 in Hg)		Chilltip contact cooling (2°C)		Chilltip contact cooling (2°C)	
On-Site Service	Yes; replaceable handpiece		Yes; replaceable handpiece		Yes; replaceable handpiece	
Connectivity	Remote connectivity for service and support					
Platform Specifications						
Electrical	100-240 VAC +/-10%, 10 A 50/60 Hz. Single-phase grounded outlet					
Weight	The console weight is: 36 kg (79.4 lbs.)					
Dimensions (WxDxH)	51 x 60 x 42 cm; (20.1 x 23.6 x 16.55 in)					

Indications for use

Treatment of benign vascular lesions, including angiomas, hemangiomas, telangiectasia, and other benign vascular lesions and leg veins. Hair removal, permanent hair reduction*, and the treatment for Pseudofolliculitis Barbae (PFB). Treatment of benign pigmented lesions. Treatment of wrinkles (for 1060nm handpieces only). Use on all skin types (Fitzpatrick skin type I-VI), including tanned skin.

* Permanent hair reduction is defined as the long-term, stable reduction in the number of hairs re-growing when measured at 6, 9 and 12 months after the completion of a treatment regime.

Trademarks

LightSheer, LightSheer QUATTRO, the LightSheer QUATTRO logo, ChillTip, Pulse Width (OptiPulse), AutoPace, Lumenis and the Lumenis logo are trademarks or registered trademarks of the Lumenis group of companies.

About Lumenis

Lumenis is a global leader in the field of minimally invasive clinical solutions for the surgical, ophthalmology and aesthetic markets, and is a world-renowned expert in developing and selling energy-based technologies, including laser, intense pulsed light (IPL) and radio frequency (RF). For over 50 years, Lumenis' ground-breaking products have redefined medical treatments and have set numerous technological and clinical gold standards.













Lumenis LTD

Yokneam Industrial Park Hakidma Street 6, P.O.B. # 240 Yokneam 2069204, ISRAEL T +972 4 959 9000

Lumenis (Germany) GmbH

Heinrich-Hertz-Str 3 D-63303 Dreieich-Dreieichenhain GERMANY T +49 (0) 6103 8335 0

EA

Dreieich Dreieichenhain, Germany T +49 6103 8335 0

F +49 6103 8335 300

Roma, Italy

T +39 06 90 75 230 F +39 06 90 75 269

Hertfordshire, UK T +44 20 8736 4110

F +44 20 8736 4119

Japan

Tokyo, Japan T +81 3 6743 8300 F +81 3 6743 8301

AMERICAS

San Jose, CA, USA T +1 408 764 3000 +1 877 586 3647 F +1 408 764 3999

ASIA PACIFIC Beijing, China

T + 86 10 5737 6677 F +86 10 5737 6767

Gurgaon, India T + 91 124 422 07 95

Kowloon, Hong Kong

T + 852 217 428 00 F + 852 272 251 51





© 2019 Lumenis Ltd., All Rights Reserved.





De Grift 20 7711 EJ, Nieuwleusen